



# Mammotome Confirm™

Core Specimen Radiography System



Instant and clear  
biopsy core  
imaging at your  
patient's side

Improve biopsy workflow with the Mammotome Confirm™ system:



### Optimized efficiency

Improve biopsy suite utilization by eliminating unnecessary sample analysis delays



### Integrated solutions

Reduce procedure time by integrating point-of-care specimen imaging directly into the biopsy suite



### Improved quality

Identify and verify even challenging microcalcifications at the time of biopsy

## System Features


- High resolution specimen images with exceptional clarity within seconds
- Expansive imaging field accommodating up to 12 specimens in a single image
- Wireless communication for easy transmission of one or more images with a single click
- Touchscreen option with a full set of image annotation tools
- Small footprint design minimizes space requirements



## Comparison of Imaging Area and Footprint



The Mammotome Confirm™ system has a smaller footprint than most core radiography systems, 30% smaller than the Faxitron CoreVision.

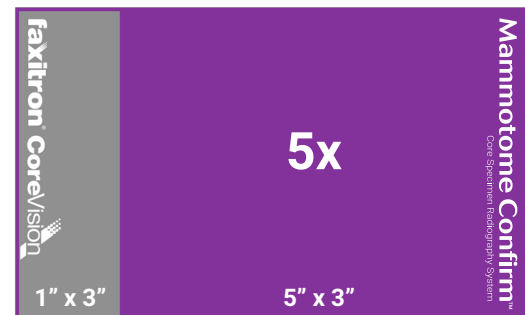

The Mammotome Confirm™ system has 2X the imaging area at 32 $\mu$  the highest resolution and 5X the imaging area at standard resolution.

Item	Specification Description
Detector Size	3" x 5" (12 cm x 7 cm)
Spatial Resolution	10.4 lp/mm, contact mode   15.6 lp/mm, mag mode
Energy Range	10-25 kV
Tube Current	up to 1.0 mA
Window Filtration	0.005" beryllium
Focal Spot	50 $\mu$ m, nominal
Power	110/220 VAC, 50/60 Hz, 500 VA Max
Wireless Connectivity	Standard
DICOM Interface	Annotate, Store, Print, and Modality Work List
Size (W x D x H)	11" x 12" x 20" (28 x 30.5 x 51 cm)
Weight Without cart	70 lbs. (32 kg)
Weight With cart	170 lbs. (77 kg)

\*cart sold separately

Contact your local Mammotome Representative to learn more

1-877-926-2666  
mammotome.com

**Mammotome**

Product may not be approved or available in your region. Please check with your local Mammotome representative.

© 2021 Devicor Medical Products, Inc. All rights reserved. Devicor, Mammotome and Confirm are registered trademarks of Devicor Medical Products, Inc., in the USA and optionally in other countries. Other logos, product and/or company names might be trademarks of their respective owners. MDM# 1807 Rev 8/21